



## Product Sheet

# pCMV<sub>lux</sub> Vector

### Intended Use

This product is intended for research use only. It is not intended for any animal or human therapeutic or diagnostic use.

### Disclaimers

While 490 BioTech makes every effort to include accurate and up to date information on this product information sheet, 490 BioTech makes no warranties or representations of its accuracy. This product is sent with the condition that you are responsible for its safe storage, handling, and use. 490 BioTech is not liable for any damages or injuries arising from the receipt and/or use of this product.

### Limited Use Label License

This product is to be used for internal, non-commercial research purposes for the sole benefit of the purchaser. It may not be used for any other purpose, including, but not limited to diagnostics or therapeutics, and may not be used in humans. This product may not be transferred or sold to third parties, resold, modified for resale, or used to manufacture or develop commercial products or to provide a service of any kind to third parties, including, without limitation, reporting the results of purchaser's activities for a fee or other form of consideration. For information on obtaining commercial research or additional rights, please contact 490 BioTech at [licensing@490biotech.com](mailto:licensing@490biotech.com).

490 BioTech  
2450 E.J. Chapman Drive  
Knoxville, TN 37996  
[www.490BioTech.com](http://www.490BioTech.com)



### Description

The pCMV<sub>lux</sub> vector encodes a fully human expression-optimized synthetic luciferase reporter gene cassette based on the bacterial luciferase (*Photobacterium luminescens*) operon under the control of a strong cytomegalovirus (CMV) promoter for constitutive autoluminescent expression across a wide range of mammalian cell types.

### This vector contains the following features:

- Human-optimized *lux* reporter cassette for expression in mammalian cells
- Inter-gene linker regions for coordinated cassette expression
- CMV promoter for high translational expression
- Neomycin resistance gene for mammalian cell selection of the plasmid
- Ampicillin resistance gene for bacterial selection of the plasmid


**Storage Buffer:** The pCMV<sub>lux</sub> vector is supplied either in 1× TE buffer (10 mM Tris-HCl, 1 mM EDTA, pH 8.0) or in a lyophilized (dried) format. If received lyophilized, reconstitute by centrifuging at 6,000 × g for 1 minute to pellet the DNA to the bottom of the tube, adding 50 µl nuclease-free water or TE buffer, vortexing for 1 minute, and incubating at room temperature for 5 minutes (If necessary, heat at 50° C for 5-10 minutes to fully dissolve the DNA into solution). The vial can then be briefly vortexed and pulse centrifuged to collect the re-suspended DNA to the bottom of the tube and the reconstituted DNA will be ready for use.

**Storage Conditions:** Store at -20° C. Avoid multiple freeze-thaw cycles and exposure to frequent temperature changes. Temperature fluctuations can greatly alter product stability.

**Usage Note:** Concentration gradients may form in frozen products and should be dispersed upon thawing. Mix well prior to use.

### pCMV<sub>lux</sub> Feature List:

CMV immediate early enhancer/promoter	30 - 780 bp
Human optimized synthetic luciferase cassette	971 - 7675 bp
Neomycin resistance region	8450 - 9814 bp
Bacterial replication origin	10039 - 10720 bp
Ampicillin resistance region	10817 - 11477 bp

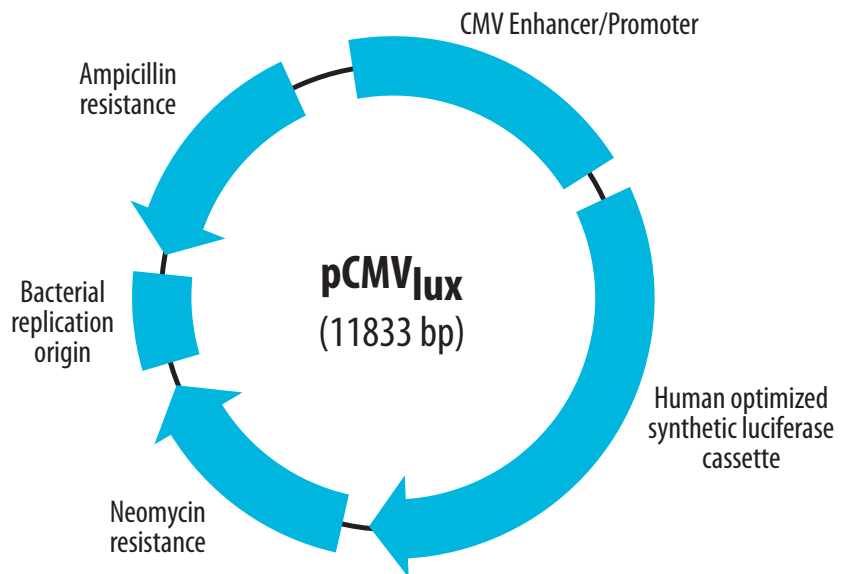

**Restrictions on Modification:** Researchers shall have no right to modify or otherwise create variations of the nucleotide sequence of the luciferase cassette genes except that researchers may: (1) create fused gene sequences provided that the coding sequence of the resulting luciferase cassette has no more than four deoxynucleotides missing at the affected terminus compared to the intact luciferase cassette sequence, and (2) insert and remove nucleic acid sequences in splicing research predicated on the inactivation or reconstruction of the luminescence of the encoded luciferase cassette. No other use or transfer of this product or derivatives is authorized without the prior express written consent of 490 BioTech.

490 BIOTECH MAKES NO REPRESENTATIONS OR WARRANTIES OF ANY KIND, EITHER EXPRESSED OR IMPLIED, INCLUDING FOR MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE WITH REGARDS TO THIS PRODUCT.

*The substrate free revolution.*